

**Product Name:** CRBN-BRD4 Proximity/Competition Assay Kit

**Product Code:** TAK-107

**FOR RESEARCH USE ONLY (RUO)**

## Shipping & Storage


**Shipping:** The product is shipped with dry ice. Upon receipt, store it immediately at the temperature recommended below.

**Stability/Storage:** Use a manual defrost freezer and avoid repeated freeze-thaw cycles. Store  $\leq -70^{\circ}\text{C}$  (stable for 24 months from date of receipt).

## Background

### Description

The CRBN-BRD4 Proximity/Competition Assay Kit provides a proximity-based screening assay for identifying compounds that bind the canonical CRBN ligand-binding site. Using dBET1 to induce a CRBN-BRD4 complex and Alpha-compatible detection beads, compounds that compete for CRBN binding reduce the proximity signal and can be readily detected in a 384-well format.

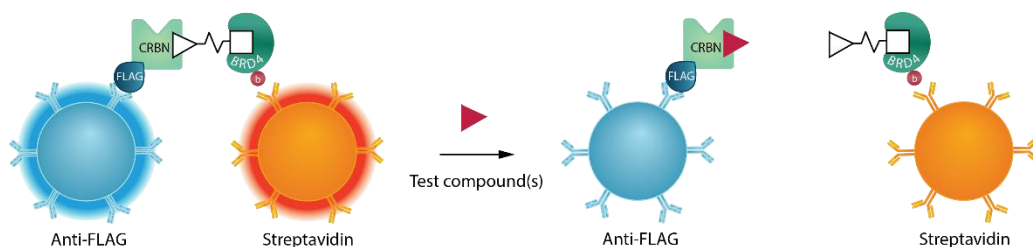

## Related Products

**TE3-015:** CRBN / DDB1

**TAP-077:** BRD4-BD1

**TBS-089:** 5X Proximity Assay Buffer