

Product Name: VHL-BRD4 Proximity/Competition Assay Kit

Product Code: TAK-108

FOR RESEARCH USE ONLY (RUO)

Shipping & Storage


Shipping: The product is shipped with dry ice. Upon receipt, store it immediately at the temperature recommended below.

Stability/Storage: Use a manual defrost freezer and avoid repeated freeze-thaw cycles. Store $\leq -70^{\circ}\text{C}$ (stable for 24 months from date of receipt).

Background

Description

The VHL-BRD4 Proximity/Competition Assay Kit provides a proximity-based screening assay for identifying compounds that bind the canonical VHL ligand-binding site. Using the VHL-recruiting degrader MZ1 to induce a VHL-BRD4 complex and Alpha-compatible detection beads, compounds that compete for VHL binding reduce the proximity signal and can be readily detected in a 384-well format.

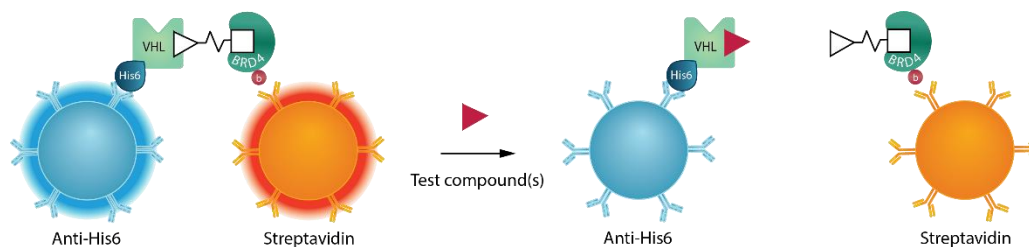

Related Products

TE3-019: VHL / ELOB / ELOC

TAP-077: BRD4-BD1

TBS-089: 10X Proximity Assay Buffer